

Name of Child: _____ Date: _____

How to Cut Earmold Tubing for Your Child's Hearing Aids

When you get new earmolds for your child's hearing aids, you will need to cut the tubing to fit correctly. Follow these steps to cut the tubing to the correct length:


- 1) Put the new earmold with the long tubing into your child's ear.
- 2) Hang the hearing aid on the ear the way it normally sits.
- 3) Line up the tubing with the earhook of the hearing aid.
- 4) Use a pen to mark the tubing just above where it connects to the notched part of the ear hook.
- 5) Remove the earmold from your child's ear and cut the tubing at the mark. **Important:** It is better to cut the tubing a little longer than needed rather than cutting it too short. You can always trim off more if needed.
- 6) Attach the trimmed earmold to the hearing aid. **Important:** Make sure the correct hearing aid is connected to the correct earmold. The right earmold should be connected to the hearing aid with a red marking. The left earmold should be connected to the hearing aid with a blue marking.



1) Put the new earmold with the long tubing into your child's ear



4) Use a pen to mark the tubing

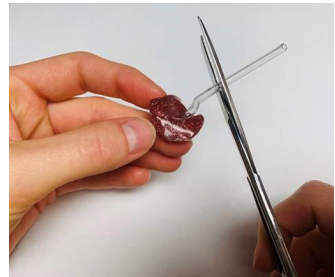

5) Cut the tubing at the mark



6) Attach the trimmed earmold to the hearing aid.

- 7) Put the hearing aid on your child's ear. The hearing aid with the red marking goes on the right ear. The hearing aid with the blue marking goes on the left ear. Check to make sure the hearing aid fits correctly:
- If the tubing is too long, it will move around or sit off the ear.
 - If the tubing is too short, it will feel too tight and may leave a red mark.
 - If you cut the tubing too short, call your child's audiologist so they can re-tube the earmold for you.



7) Put the hearing aid on your child's ear

Now that you have read this:

- Show your child's audiologist how you will trim the earmold tubing for the hearing aids. (Check when done.)



If you have any questions or concerns,
call the Phoenix Children's Audiology Department at
602-933-0905

If you want to know more about child health and illness,
visit our library at The Emily Center at Phoenix Children's Hospital
1919 East Thomas Road
Phoenix, AZ 85016
602-933-1400
866-933-6459
www.phoenixchildrens.org
www.theemilycenter.org
Facebook: [facebook.com/theemilycenter](https://www.facebook.com/theemilycenter)
Twitter: @emilycenter

Disclaimer

The information provided at this site is intended to be general information, and is provided for educational purposes only. It is not intended to take the place of examination, treatment, or consultation with a physician. Phoenix Children's Hospital urges you to contact your physician with any questions you may have about a medical condition.

November 12, 2020 • In family review
#3121 • Written by Wendy Steuerwald, Au.D, MS-ELOC

How to Cut Earmold Tubing for Your Child's Hearing Aids

Family Review of this Handout

Families: Please let us know what you think of this handout.

Would you say this handout is easy to read? Yes No
Please circle the parts of the handout that were hard to understand.

Would you say this handout is interesting to read? Yes No
— Why or why not?

Would you do anything differently after reading this handout? Yes No
— If yes, what?

After reading this handout, do you have any questions about the subject? Yes No
— If yes, what?

Is there anything you do not like about the drawings or pictures? Yes No
— If yes, what?

What changes would you make in this handout to make it better or easier to understand?

Please return your review of this handout to your nurse or doctor or send it to the address below. You can also scan it and email it to emilycenter@phoenixchildrens.com.

The Emily Center
Health Education Specialist
Phoenix Children's Hospital
1919 East Thomas Road
Phoenix, AZ 85016-7710